


Analysis of Glimepiride

(Under the Condition of the Japanese Pharmacopoeia)

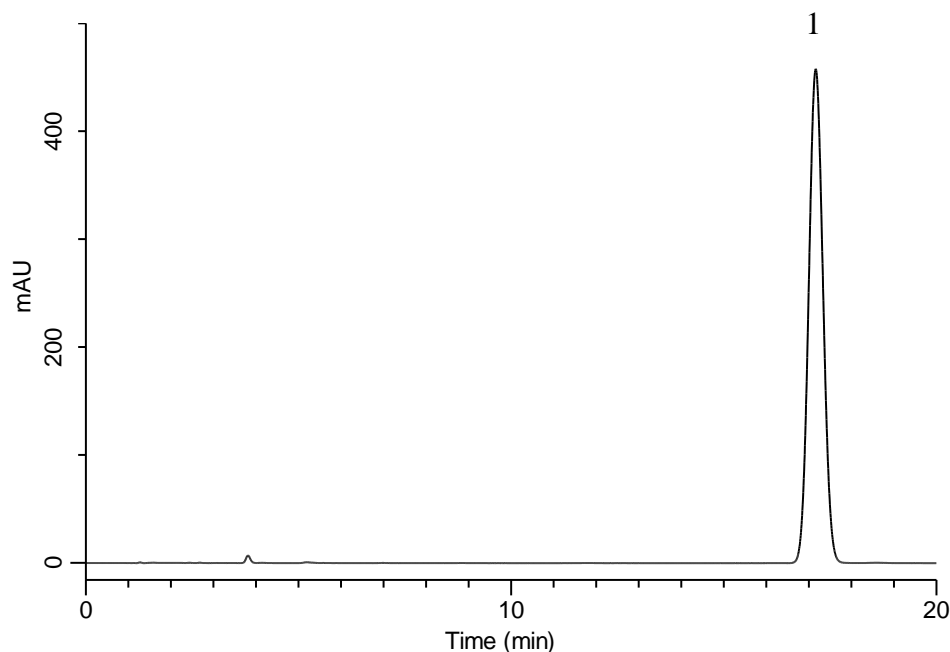

Conditions

System : GL-7400 HPLC system
Column : Inertsil ODS-3
 (4 μ m, 250 x 4.0 mm I.D.)
Column Cat. No. : 5020-04636
Eluent : A) CH₃CN
 B) 1.0 g/L NaH₂PO₄·2H₂O (pH 2.5, H₃PO₄)
 A/B = 50/50, v/v
Flow Rate : 1.15 mL/min
Col. Temp. : 25 °C
Detection : UV 228 nm (GL-7452A PDA Detector)
Injection Vol. : 20 μ L

Analyte:

1. Glimepiride 200 mg/L

Theoretical plates : 12,254 (\geq 9,000)
 Tailing factor : 1.05 (\leq 1.5)
 RSD of the
 peak area (%) (n=6) : 0.70 (\leq 1.0)