


## Analysis of Steroidal anti-inflammatory drug

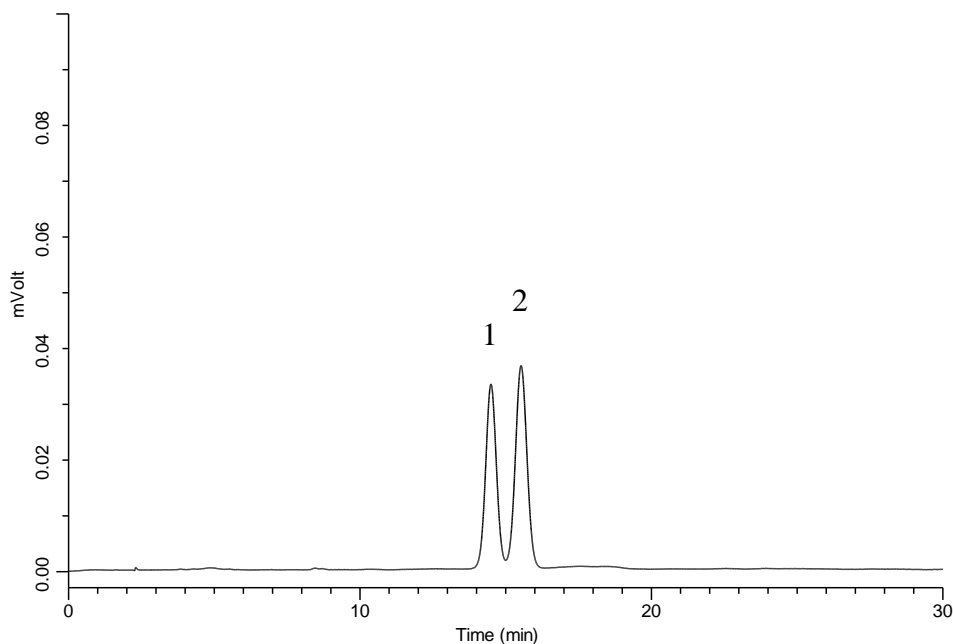

### Conditions

**Column** : InertSustain Phenyl  
(5  $\mu$  m, 150 x 4.6 mm I.D.)

**Column Cat. No.** : 5020-16327

**Eluent** : A) CH<sub>3</sub>OH  
B) H<sub>2</sub>O  
A/B = 40/60, v/v

**Flow Rate** : 0.8 mL/min

**Col. Temp.** : 40 °C

**Detection** : UV 220 nm

**Injection Vol.** : 5  $\mu$  L

**Sample** : Standard

### Analyte:

1. Hydrocortisone
2. Prednisolone