


InertSearch™ for LC

Inertsil® Applications

Analysis of Tulathromycin

Data No. LB271-0776

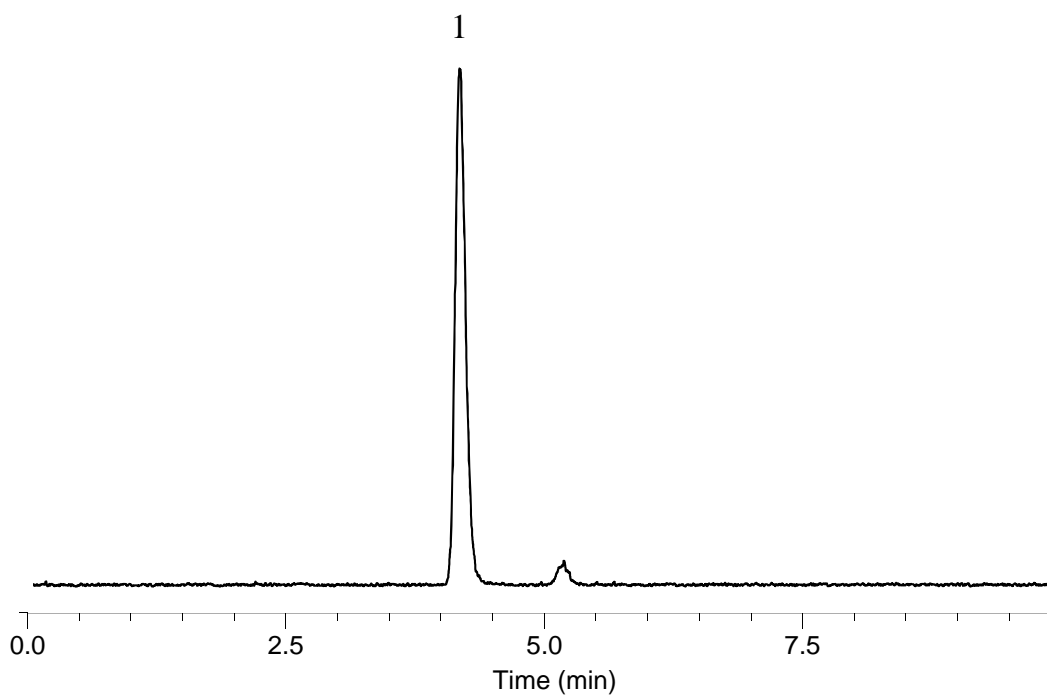

Conditions

Column : InertSustainSwift C18
(3 μ m, 150 x 2.1 mm I.D.)
Column Cat. No. : 5020-88128
Eluent : A) CH₃CN
B) 0.05 % TFA in H₂O
A/B = 25/75, v/v
Flow Rate : 0.2 mL/min
Col. Temp. : 40 °C
Detection : LC/MS/MS (ESI, Positive, MRM)
Injection Vol. : 1 μ L
Sample : Standard

Analyte:

1. Tulathromycin 10 μ g/L