


# InertSearch™ for LC

Inertsil® Applications

## Analysis of *o*-Phthalaldehyde

(Directed by Pre-Column Method with DNPH)

Data No. LB326-0659

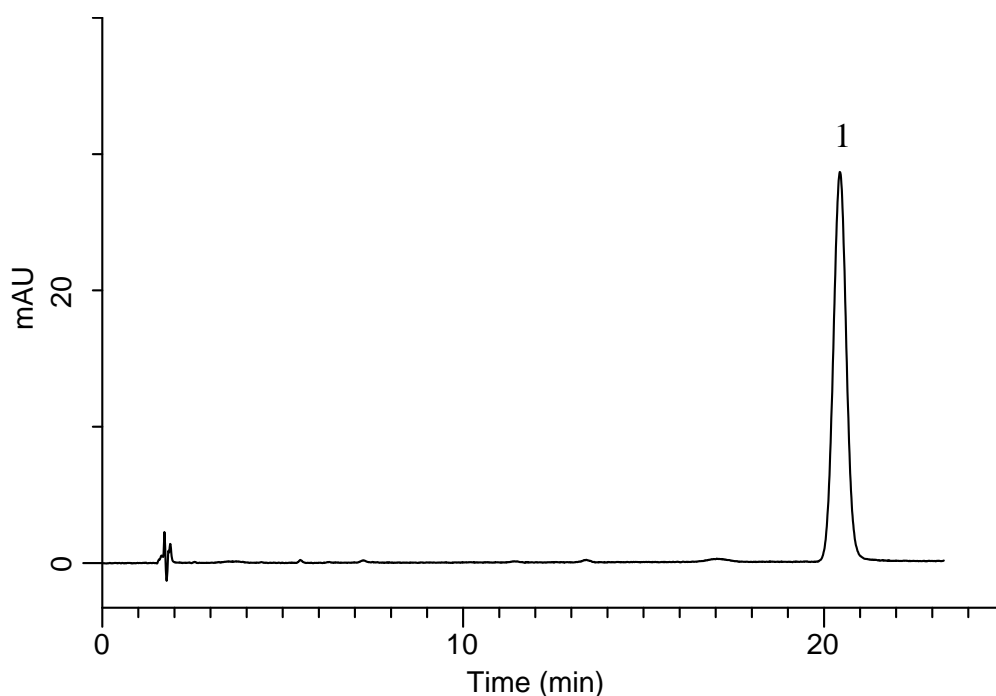

### Conditions

**System** : GL7700 HPLC system  
**Column** : InertSustainSwift C18  
(5  $\mu$  m, 150 x 4.6 mm I.D.)  
**Column Cat. No.** : 5020-88026  
**Eluent** : A) CH<sub>3</sub>CN  
B) H<sub>2</sub>O  
A/B = 60/40, v/v  
**Flow Rate** : 1.0 mL/min  
**Col. Temp.** : 40 °C  
**Detection** : UV 360 nm (PD7752 PDA Detector)  
**Injection Vol.** : 5  $\mu$  L  
**Sample** : Standard

### Analyte:

1. DNPH-*o*-Phthalaldehyde  
(10  $\mu$  g/mL in CH<sub>3</sub>CN/DMSO = 70/30, v/v)