


# InertSearch™ for LC

Inertsil® Applications

## Analysis of Vitamin D

Data No. LB467-0812

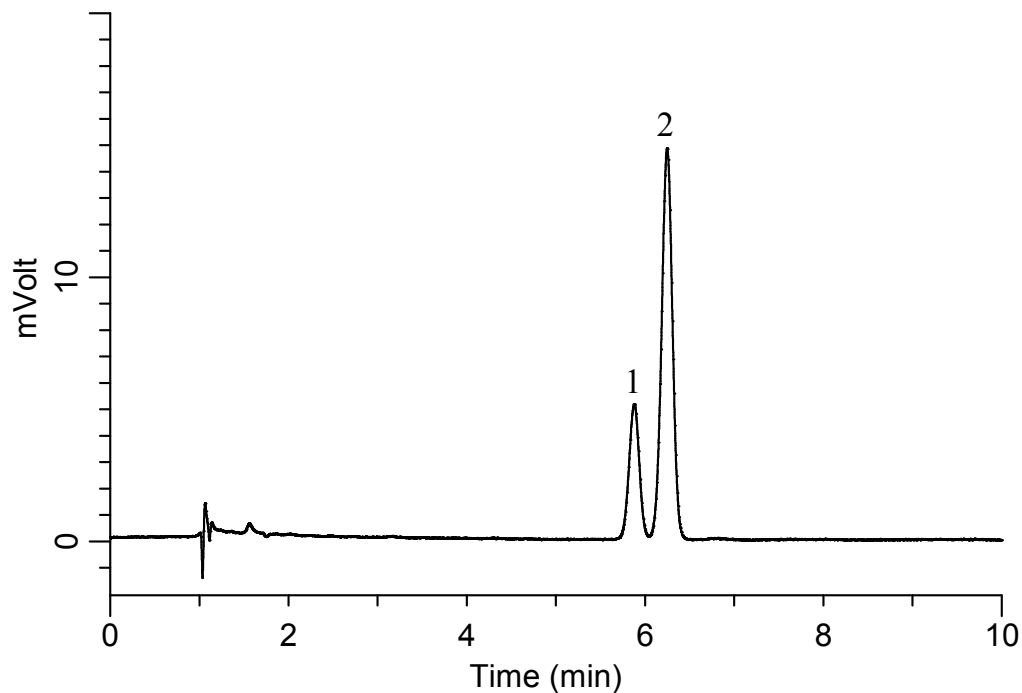

### Conditions

**Column** : Inertsil ODS-HL  
(3  $\mu$  m, 150 x 2.1 mm I.D.)  
**Column Cat. No.** : 5020-87231  
**Eluent** : CH<sub>3</sub>OH  
**Flow rate** : 0.3 mL/min  
**Col. Temp.** : 25 °C  
**Detection** : UV 265 nm  
**Injection Vol.** : 5  $\mu$  L  
**Sample** : Standard

### Analyte:

1. Vitamin D2  
2. Vitamin D3  
(5 mg/L each)