


## Analysis of Eugenol

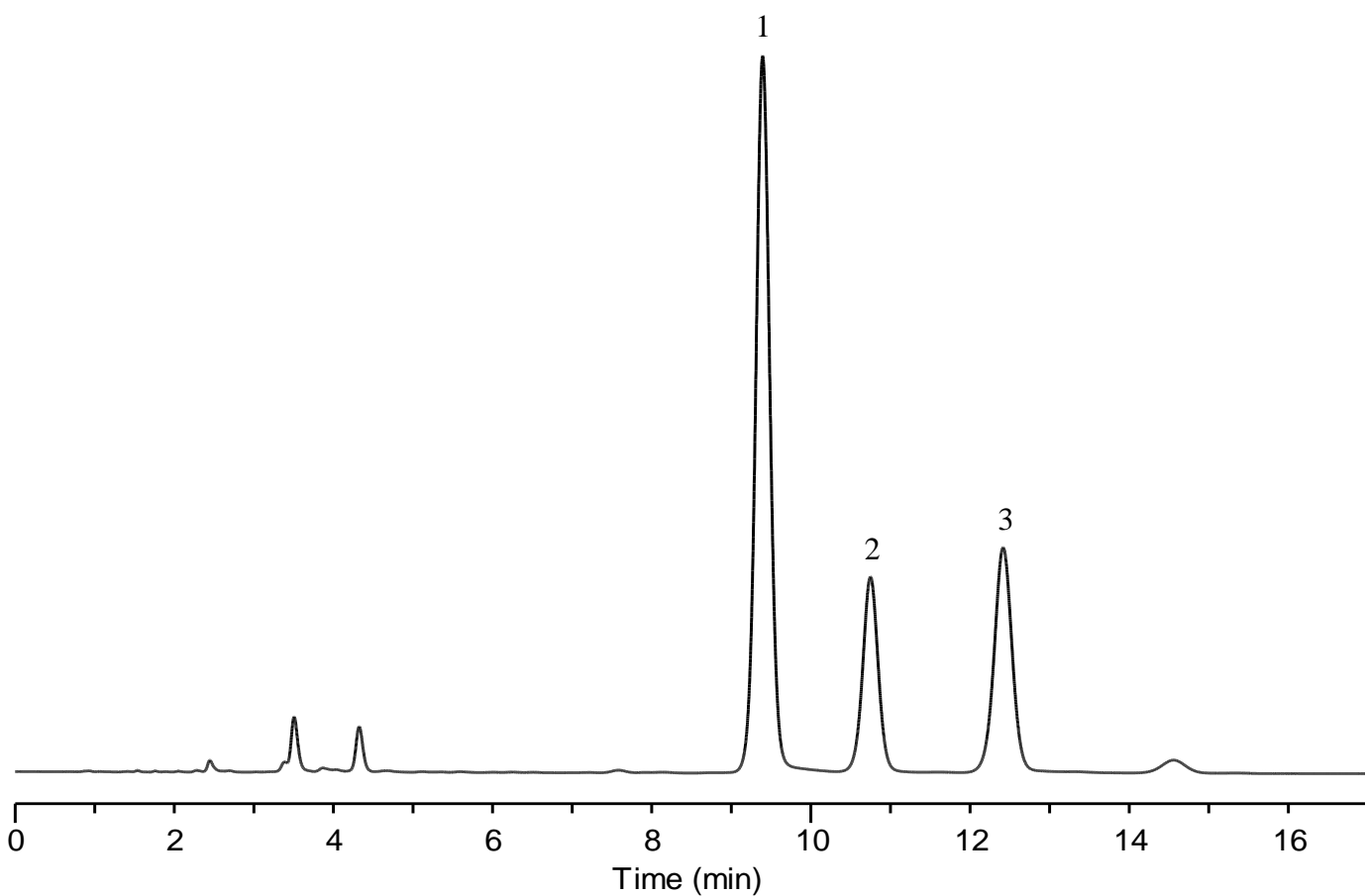

### Conditions

**System** : GL7700 HPLC system  
**Column** : InertSustain PFP HP (3  $\mu$  m, 150 x 2.1 mm I.D.)  
**Column Cat. No.** : 5020-87921  
**Eluent** : A) CH<sub>3</sub>OH  
           B) H<sub>2</sub>O  
           A/B = 40/60, v/v  
**Flow rate** : 0.2 mL/min  
**Col. Temp.** : 40 °C  
**Detection** : UV 210 nm (UV7750 with semi-micro flow cell)  
**Injection Vol.** : 2.5  $\mu$  L  
**Sample** : Standard

### Analyte:

1. Eugenol	100 mg/L
2. <i>cis</i> -Isoeugenol	100 mg/L
3. <i>trans</i> -Isoeugenol	100 mg/L