


Analysis of Sertraline

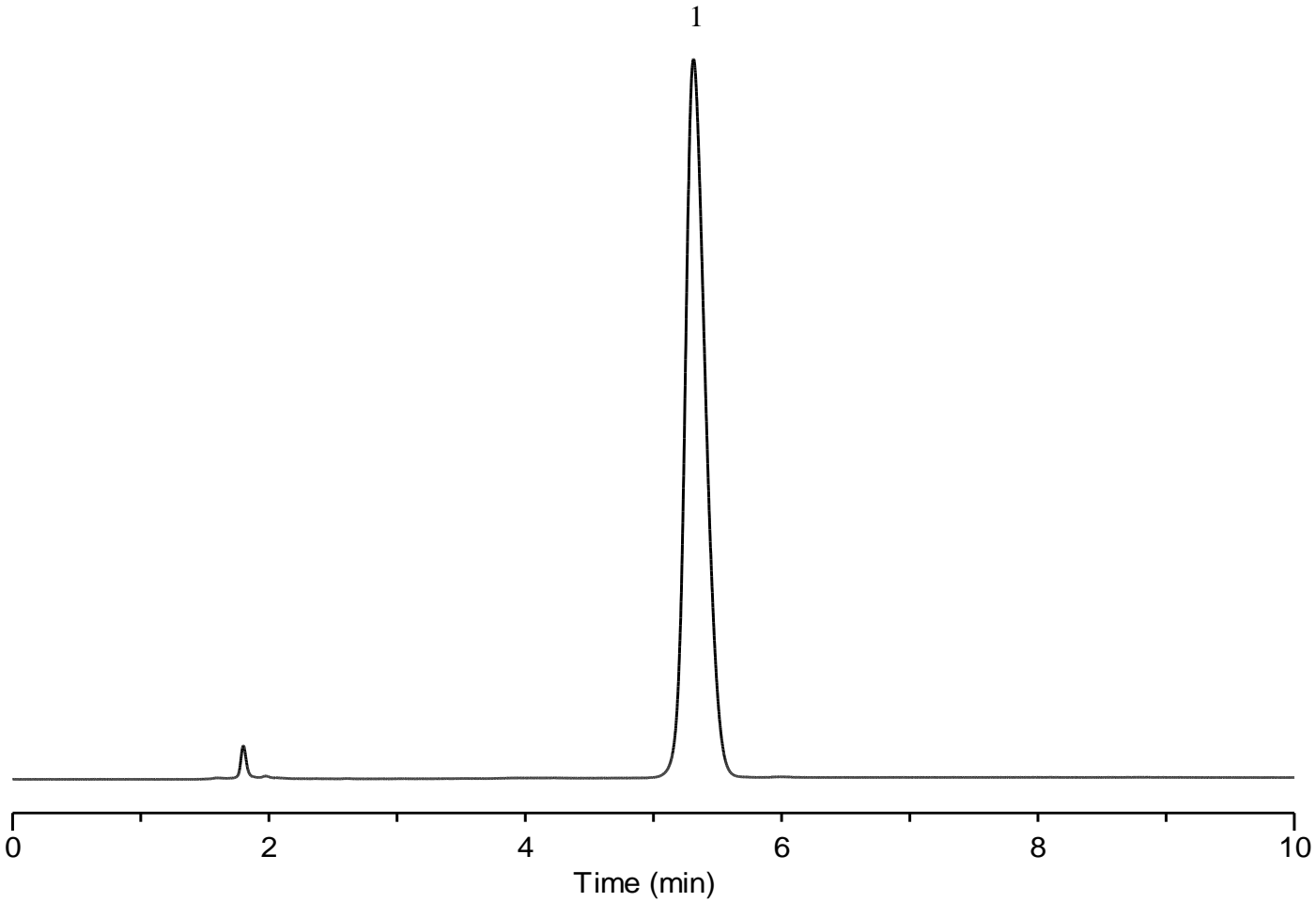

Conditions

Column : InertSustain Cyano (5 μ m, 250 x 4.6 mm I.D.)
Column Cat. No. : 5020-89277
Eluent : A) CH₃OH
 : B) 0.1% H₃PO₄ in H₂O
 : A/B = 50/50, v/v
Flow rate : 1.5 mL/min
Col. Temp. : 30 °C
Detection : UV 210 nm
Injection Vol. : 10 μ L
Sample : Standard

Analyte:

1. Sertraline 50 mg/L
 Symmetry factor (1) : 1.24