


InertSearch for LC

Inertsil Applications

Analysis of Ribonuclease A

Data No. LB667-0812

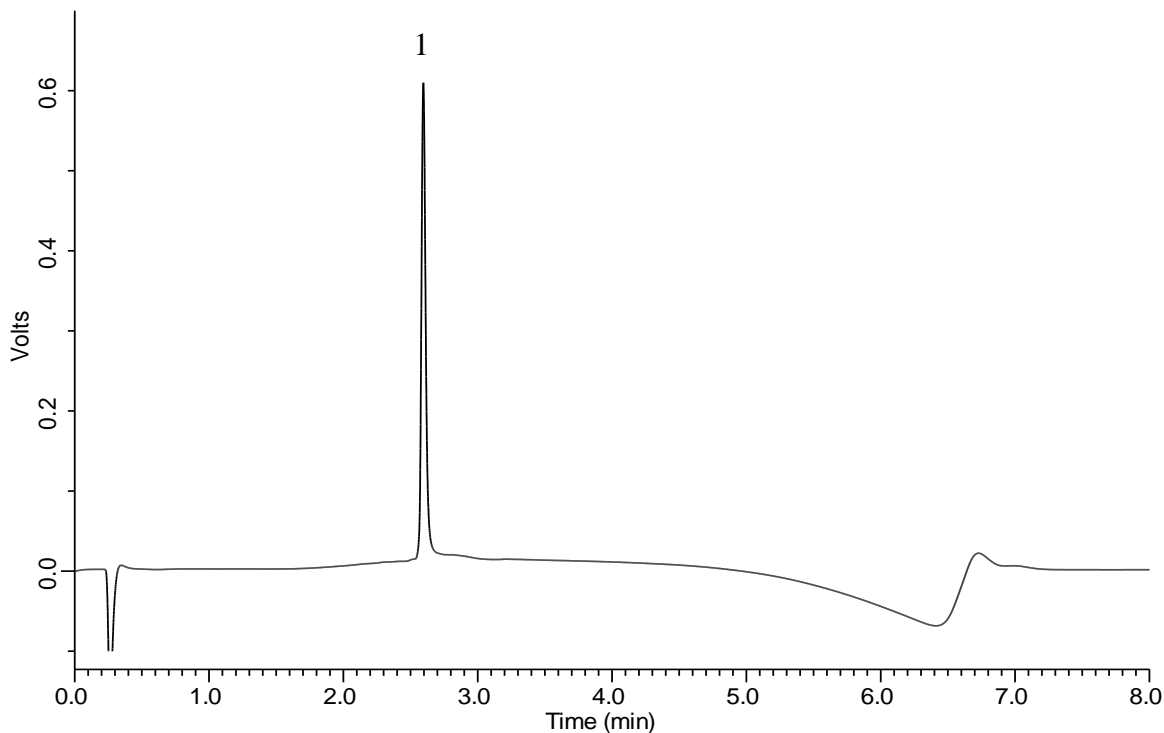

Conditions

Column Cat. No. : 5020-10818
Column : MonoSelect RP-mAb
(20 x 2.1 mm I.D.)
Eluent : A) 0.1% TFA in CH₃CN
B) 0.1% TFA in H₂O
A/B = 5/95-5 min-90/10-0.1 min-5/95-3 min
Flow Rate : 0.3 mL/min
Col. Temp. : 80 °C
Detection : UV 210 nm
Injection Vol. : 5 μL
Sample : Standard

Analyte:

1. Ribonuclease A 0.1 mg/mL