


InertSearch for LC

Inertsil Applications

Analysis of Zebularine and Gemcitabine

Data No. LB677-0919

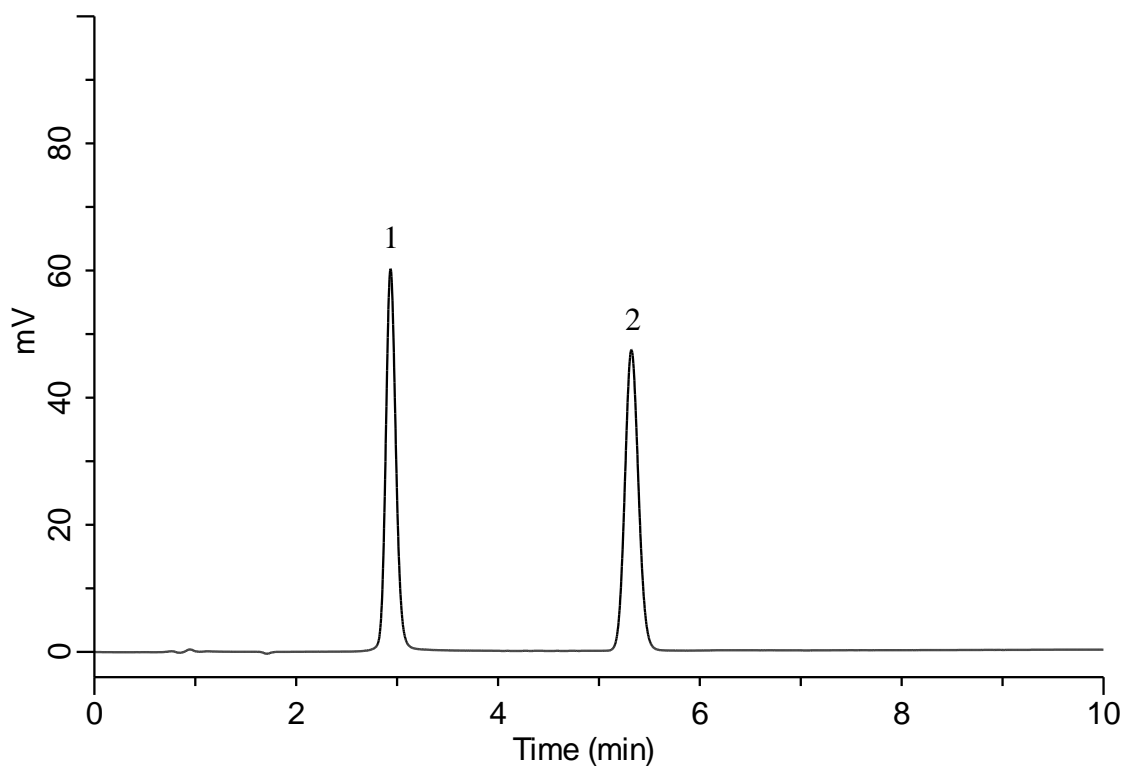

Conditions

System : Primaide
Column : InertSustain C18 HP
(3 μ m , 100 x 3.0 mm I.D.)
Column Cat. No. : 5020-14424
Eluent : 10 mM Na₂HPO₄ in H₂O (pH 7.0, H₃PO₄)
Flow Rate : 0.6 mL/min
Col. Temp. : 40 °C
Detection : UV 205 nm (Primaide 1410 UV Detector)
Injection Vol. : 10 μ L
Sample : Standard

Analyte:

1. Zebularine 10 mg/L
2. Gemcitabine 10 mg/L