


Analysis of Toluene and Toluene d8

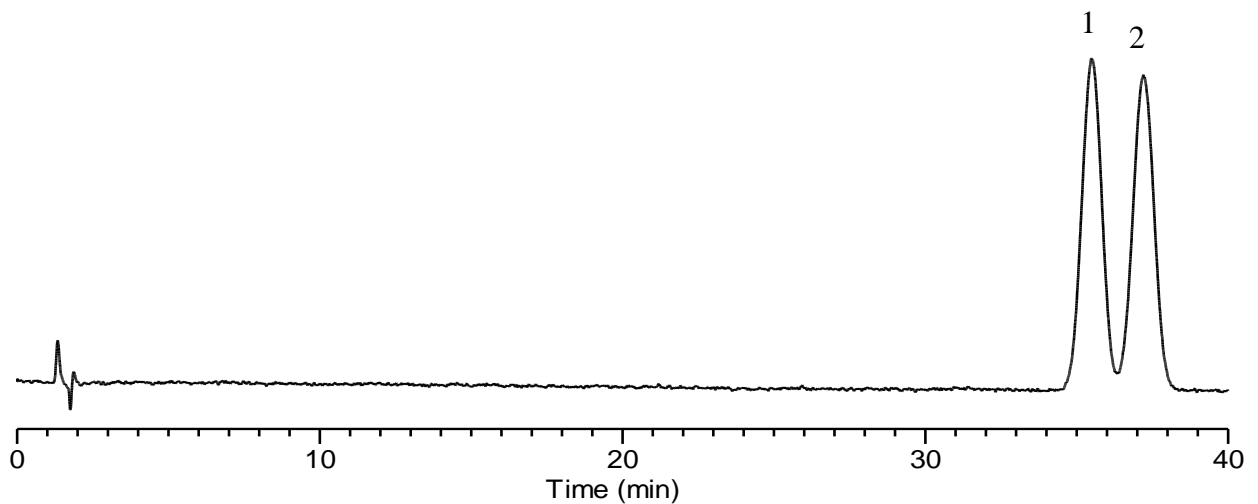

Conditions

System : Primaide HPLC system
Column : InertSustain C18 (5 μ m, 250 x 4.6 mm I.D.)
Column Cat. No. : 5020-07346
Eluent : A) CH₃CN
 : B) H₂O
 : A/B = 35/65, v/v
Flow rate : 1.5 mL/min
Col. Temp. : 40 °C
Detection : UV 254 nm
Injection Vol. : 5 μ L
Sample : Standard

Analyte:

1. Toluene-d8
2. Toluene
Conc. 100 mg/L each