


InertSearch for LC

Inertsil Applications

Analysis of Short Peptides

Data No. LB731-0806

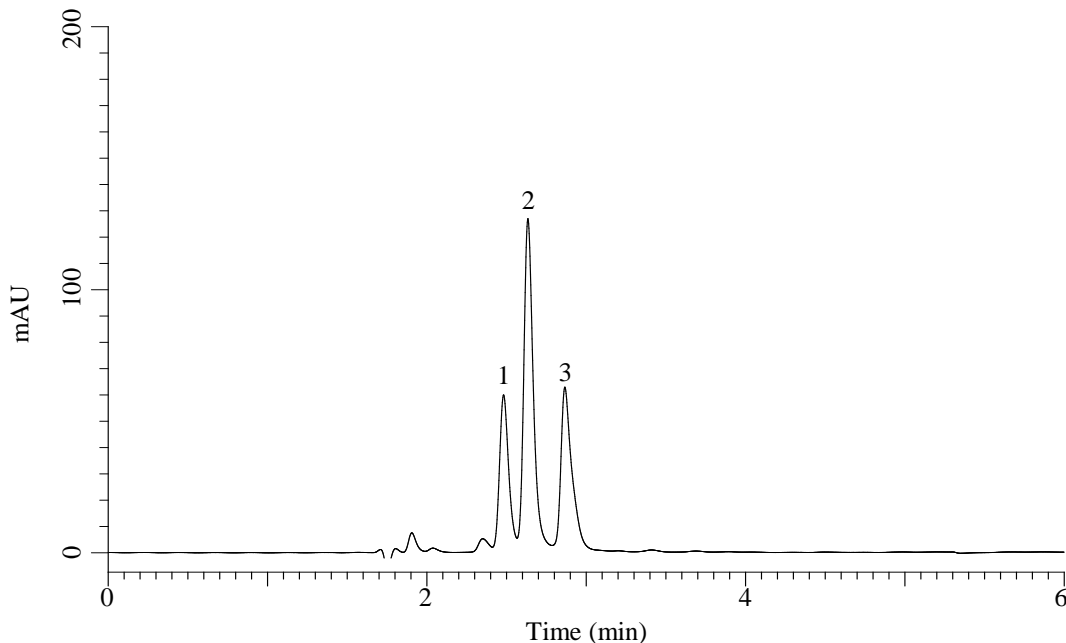

Conditions

Column : InertSustain AX-C18 (GL Sciences Inc.)
(3 μ m, 150 x 2.1 mm I.D.)
Column Cat. No. : 5020-91039
Eluent : 10 mM HCOONH₄ in H₂O
Flow Rate : 0.2 mL/min
Col. Temp. : 40 °C
Detection : UV 210 nm
Injection Vol. : 1 μ L
Sample : Standard

Analyte:

1. Gly-Arg-Gly-Asp-Ser
2. Arg-Gly-Asp
3. Arg-Gly-Asp-Ser
(250 mg/L each)