


InertSearch for LC

Inertsil Applications

Analysis of Nucleotide sugars

Data No. LB774-0806

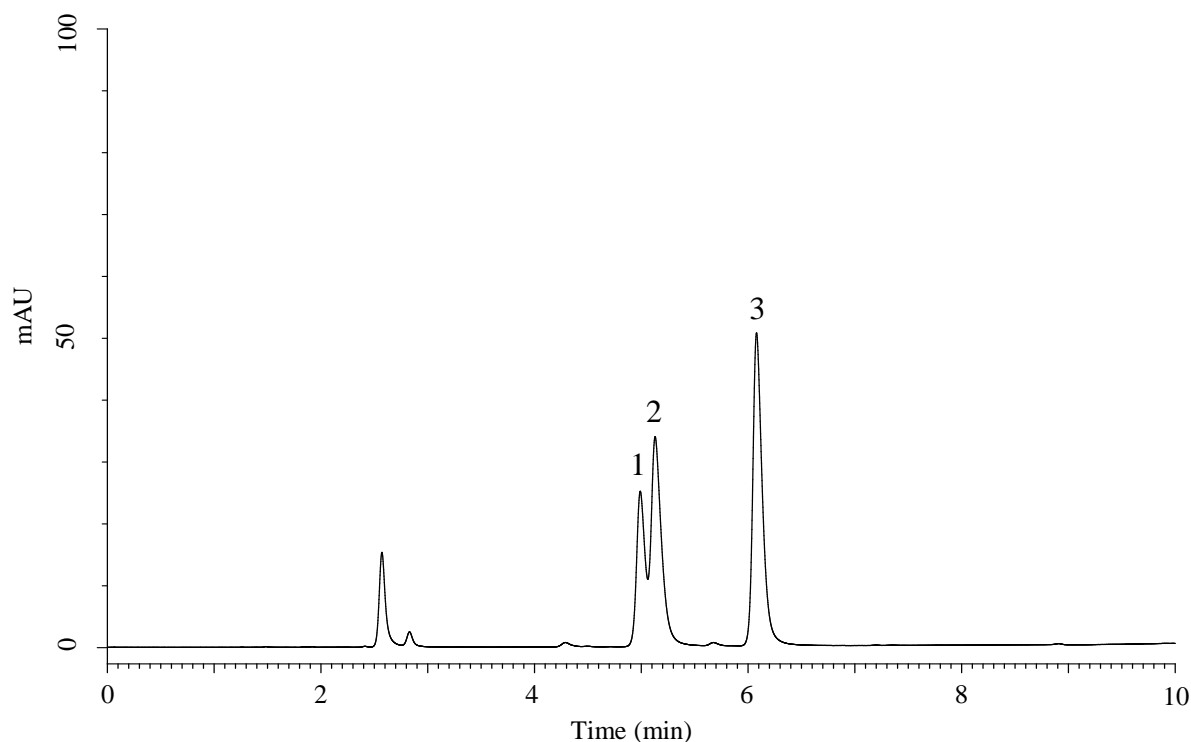

Conditions

Column : InertSustain AX-C18 (GL Sciences Inc.)
(3 μ m, 150 x 2.1 mm I.D., Metal-free hardware)

Column Cat. No. : 5020-91059

Eluent : A) 10 mM HCOONH₄ in
(CH₃CN/H₂O = 90/10, v/v)
B) 10 mM HCOONH₄ in H₂O

Time (min)	A (vol%)	B (vol%)
0.0	2	98
10.0	20	80

Flow Rate : 0.3 mL/min

Col. Temp. : 40 °C

Detection : UV 254 nm

Injection Vol. : 1 μ L

Sample : Standard

Analyte:

1. Uridine diphosphate glucose
(UDP-Glu)

2. Uridine diphosphate galactose
(UDP-Gal)

3. Guanosine diphosphate mannose
(GDP-Man)

(100 mg/L each)