


# InertSearch for LC

Inertsil Applications

## Analysis of ATP and its related compounds

Data No. LB780-0806

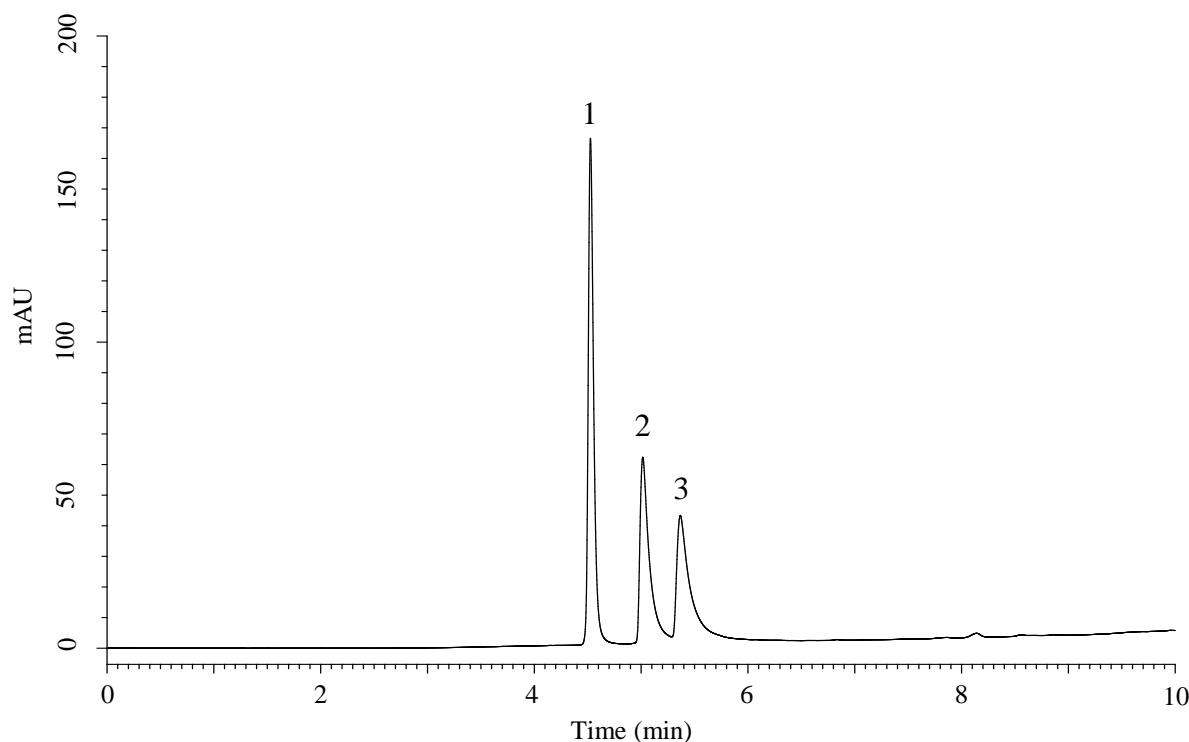

### Conditions

**Column** : InertSustain AX-C18 (GL Sciences Inc.)  
(3  $\mu$  m, 150 x 2.1 mm I.D., Metal-free hardware)

**Column Cat. No.** : 5020-91059

**Eluent** : A) 200 mM HCOONH<sub>4</sub> in CH<sub>3</sub>OH  
B) 10 mM HCOONH<sub>4</sub> in H<sub>2</sub>O

Time (min)	A (vol%)	B (vol%)
0.0	2	98
10.0	80	20

**Flow Rate** : 0.3 mL/min

**Col. Temp.** : 40 °C

**Detection** : UV 254 nm

**Injection Vol.** : 1  $\mu$  L

**Sample** : Standard

### Analyte:

1. Adenosine monophosphate (AMP)

2. Adenosine diphosphate (ADP)

3. Adenosine triphosphate (ATP)  
(100 mg/L each)