


Analysis of Vitamin D

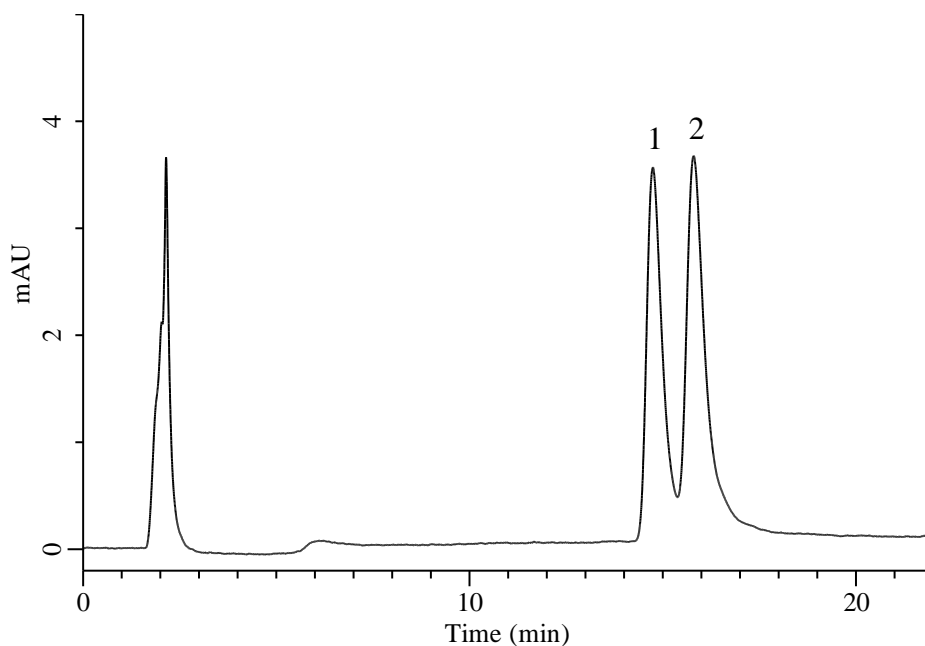

Conditions

System : Chromaster HPLC system (HITACHI)
Column : InertSustain C30 (GL Sciences Inc.)
 (HP 3 μ m, 250 x 3.0 mm I.D.)
Column Cat. No. : 5020-17181
Eluent : CH₃CN
Flow Rate : 0.6 mL/min
Col. Temp. : 25 °C
Detection : UV 265 nm
Injection Vol. : 5 μ L
Sample : Standard

Analyte:

1. Vitamin D2 5 mg/L
 2. Vitamin D3 5 mg/L