


Analysis of 5-Fluorouracil and 5-Fluorocytosine

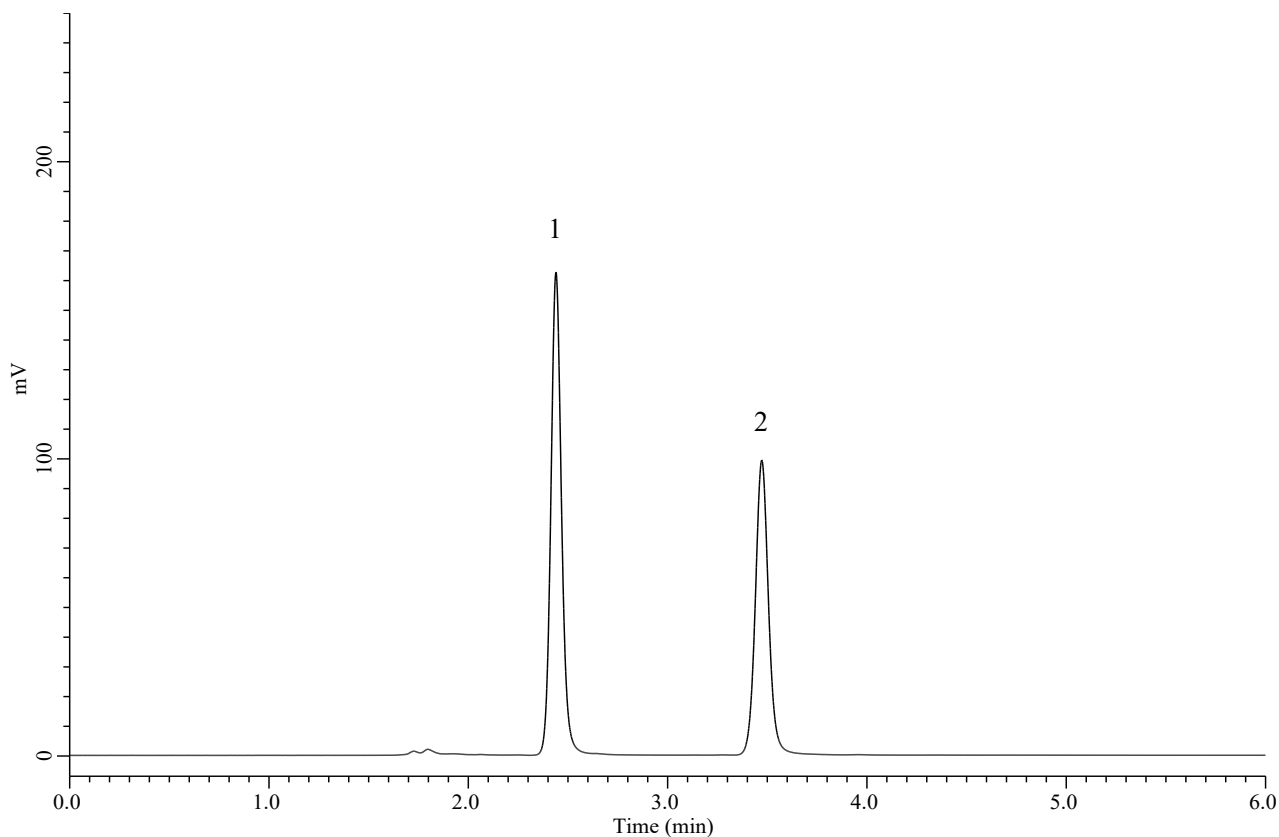

Conditions

Column : Inertsil Hybrid-Z HILIC
 (3.5 μ m, 150 x 2.1 mm I.D)
Column Cat. No. : 5020-76196
Eluent : A) CH₃CN
 B) 10 mM HCOONH₄ in H₂O
 A/B = 90/10, v/v
Flow Rate : 0.2 mL/min
Col. Temp. : 40 °C
Detection : UV 254 nm
Injection Vol. : 0.5 μ L

Analyte

1. 5-Fluorouracil
 2. 5-Fluorocytosine
 (100 mg/L of each)