


Analysis of Diquafosol

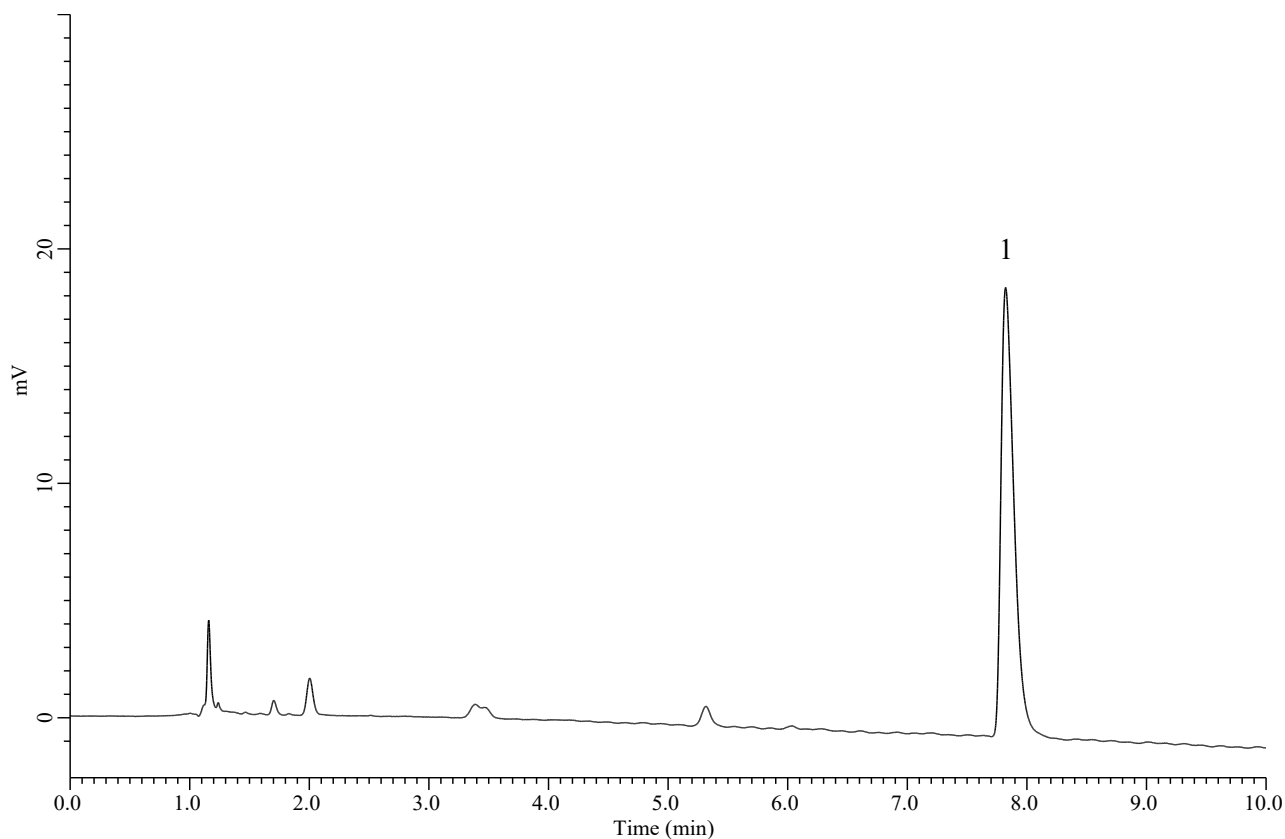

Conditions

Column : Inertsil Hybrid-Z HILIC
(3.5 μm , 150 $\times$ 2.1 mm I.D, Metal free hardware)
Column Cat. No. : 5020-76210
Eluent : A) CH_3CN
B) 10 mM HCOONH_4 in H_2O

Analyte

1. Diquafosol 50 mg/L

Time (min)	A (Vol. %)	B (Vol. %)
0	80	20
10	50	50

Flow Rate : 0.3 mL/min
Col. Temp. : 40 $^\circ\text{C}$
Detection : UV 254 nm
Injection Vol. : 1 μL