


Analysis of Erythorbic Acid and Ascorbic Acid

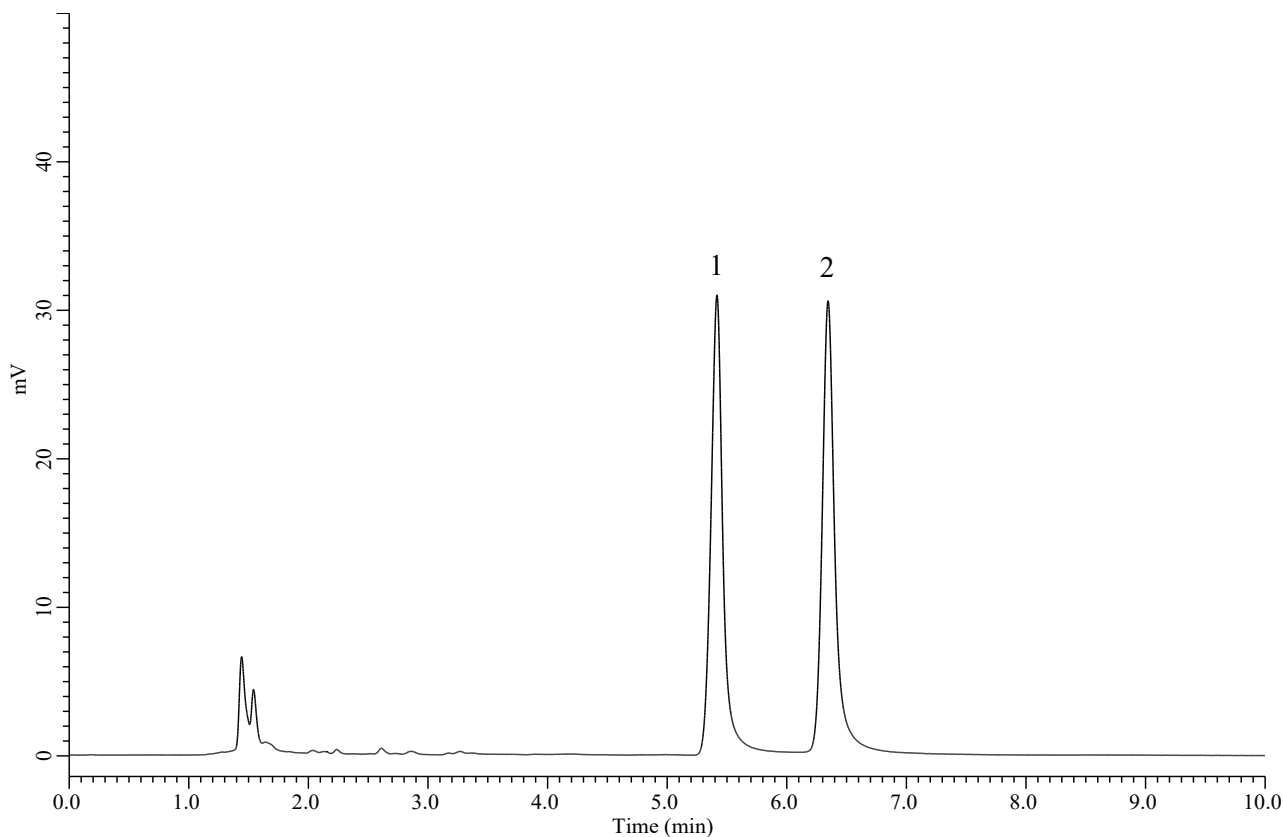

Conditions

Column : Inertsil Hybrid-Z HILIC
 (1.7 μm , 150 $\times$ 2.1 mm I.D)
Column Cat. No. : 5020-76218
Eluent : A) CH_3CN
 B) 10 mM HCOONH_4 in H_2O
 A/B = 80/20, v/v
Flow Rate : 0.2 mL/min
Col. Temp. : 40 $^\circ\text{C}$
Detection : UV 254 nm
Injection Vol. : 1 μL

Analyte

1. Erythorbic acid
 2. Ascorbic acid
 (50 mg/L of each)