


Analysis of Imidazole Dipeptide

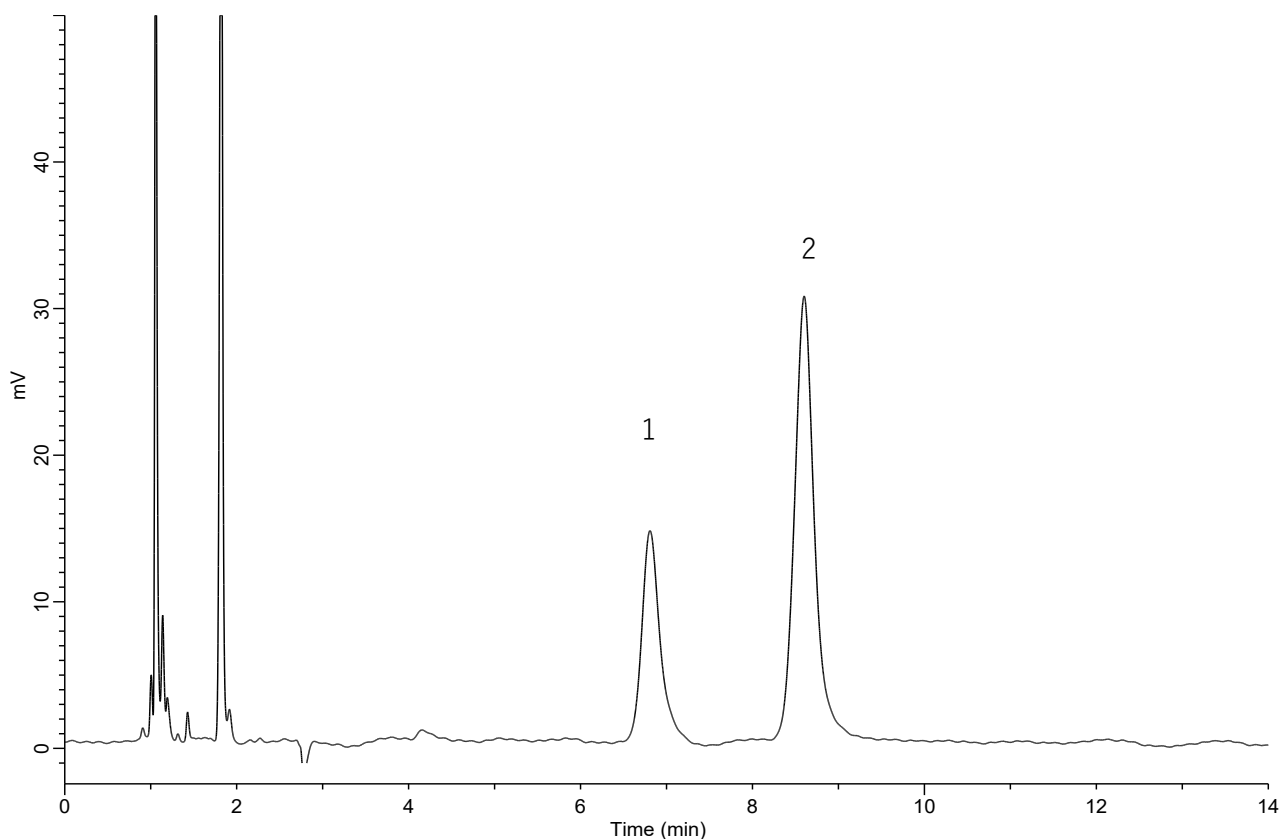

Conditions

Column : Inertsil Hybrid-Z HILIC
(3.5 μ m, 150 $\times$ 2.1 mm I.D)
Column Cat. No. : 5020-76196
Eluent : A) CH₃CN
Flow Rate : B) 0.1% HCOONH₄ in H₂O
Col. Temp. : 0.2 mL/min
Detection : 40 °C
Injection Vol. : UV 220 nm
1 μ L

Analyte

1. Carnosine
2. Anserine
(100 mg/L of each)