


Analysis of Diquafosol

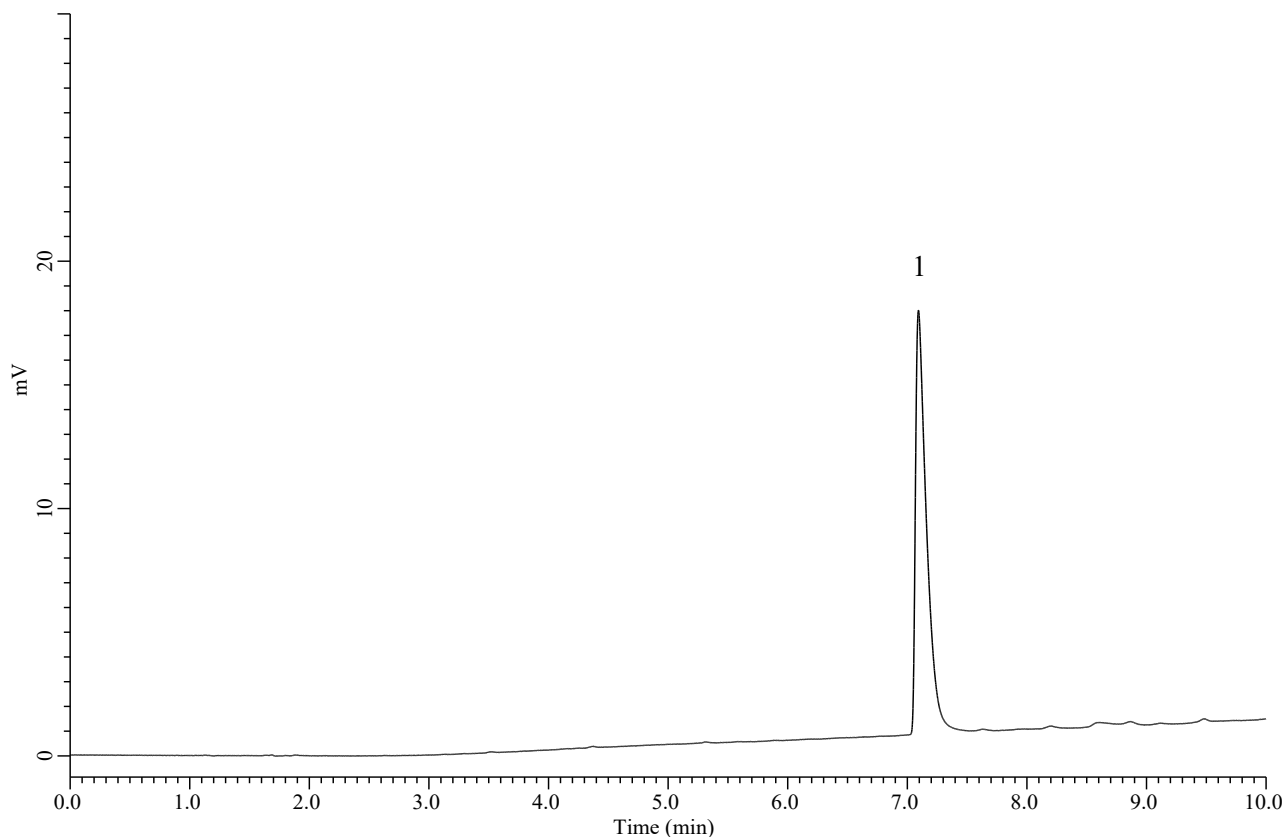

Conditions

Column : InertSustain AX-C18
(3 μm , 150 $\times$ 2.1 mm I.D, Metal free hardware)

Column Cat. No. : 5020-91059

Eluent : A) 100 mM HCOONH_4 in (H_2O : CH_3CN = 50 /50, v/v)
B) 10 mM HCOONH_4 in H_2O

Analyte

1. Diquafosol 100 mg/L

Time (min)	A (Vol. %)	B (Vol. %)
0	2	98
10	50	50

Flow Rate : 0.3 mL/min

Col. Temp. : 40 $^{\circ}\text{C}$

Detection : UV 254 nm

Injection Vol. : 0.5 μL