

# InertSearch for GC


Micro GC Applications

## CO<sub>2</sub> in H<sub>2</sub>

Data No. MA016-0000

### Specification

One GC Channel

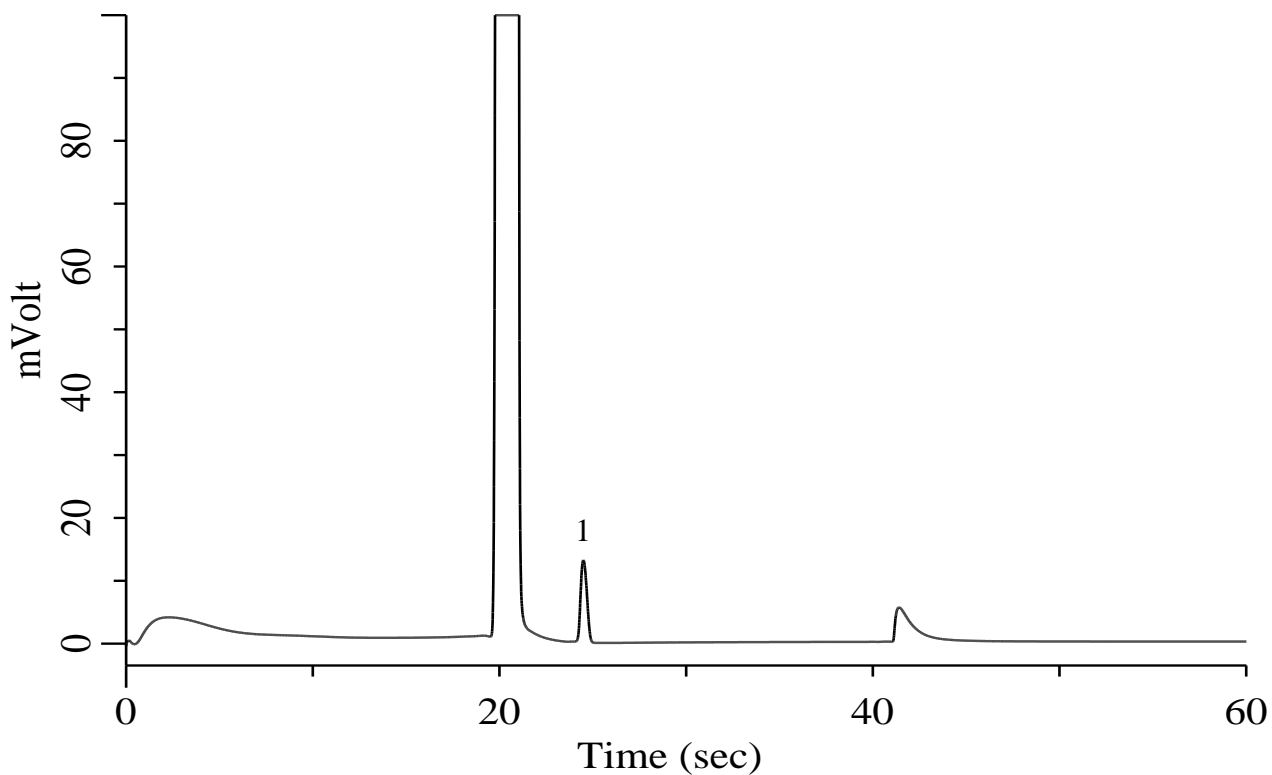

### Conditions

System : Micro GC  
Column : PoraPLOT Q 10 m  
Column Temp. : 80 °C  
Carrier Gas : Ar 190 kPa  
Inject Time : 40 msec

### Analyte:

1. Carbon dioxide (CO<sub>2</sub>) 1 % (v/v)  
Balance Gas H<sub>2</sub>