


Hydrocarbons C5-C8

Specification

One GC Channel

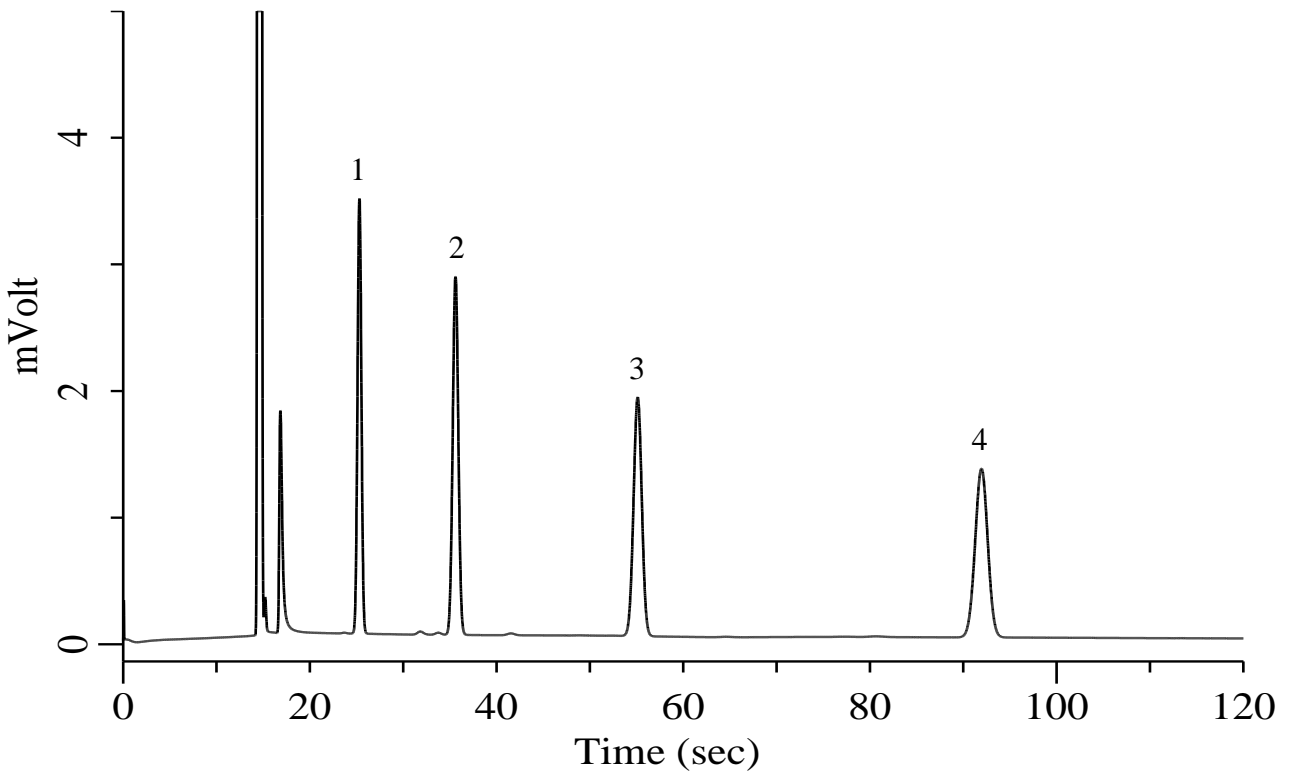

Conditions

System : Micro GC
 Column : CP-Sil 5CB 8 m
 Column Temp. : 100 °C
 Carrier Gas : He 200 kPa
 Inject Time : 40 msec

Analyte:

1. *n*-Pentane (C₅H₁₂) 1000 ppm (v/v)
 2. *n*-Hexane (C₆H₁₄) 1000 ppm (v/v)
 3. *n*-Heptane (C₇H₁₆) 1000 ppm (v/v)
 4. *n*-Octane (C₈H₁₈) 1000 ppm (v/v)

Balance Gas He