

InertSearch for GC


Micro GC Applications

Inorganic gases (N₂, Kr, Xe)

Data No. MA033-0984

Specification

One GC Channel

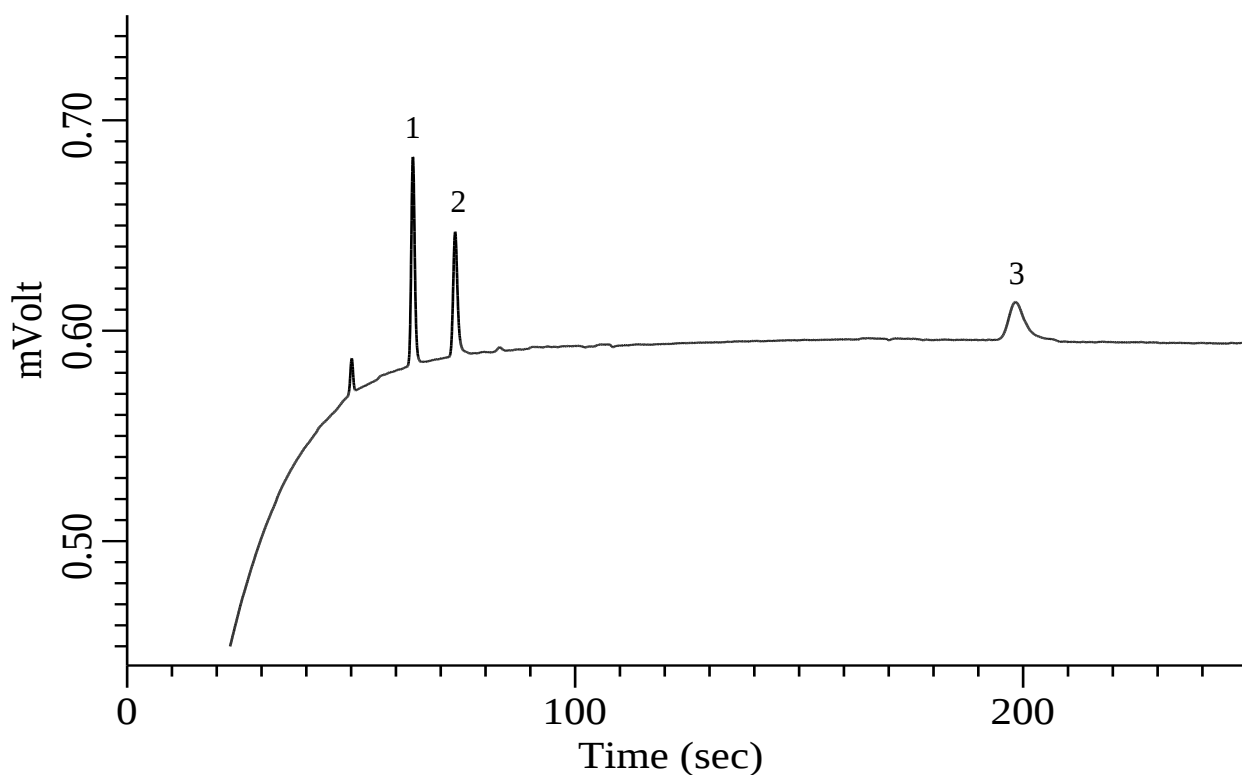

Conditions

System : Micro GC
Column : Molsieve 5A 10 m
Column Temp. : 100 °C
Carrier Gas : He 170 kPa
Inject Time : 40 msec

Analyte:

1. Nitrogen (N₂) 100 ppm (v/v)
2. Krypton (Kr) 100 ppm (v/v)
3. Xenon (Xe) 100 ppm (v/v)

Balance Gas He