

InertSearch for GC


Micro GC Applications

Organic Solvents, Hydrocarbons

Data No. MA038-0000

Specification

One GC Channel

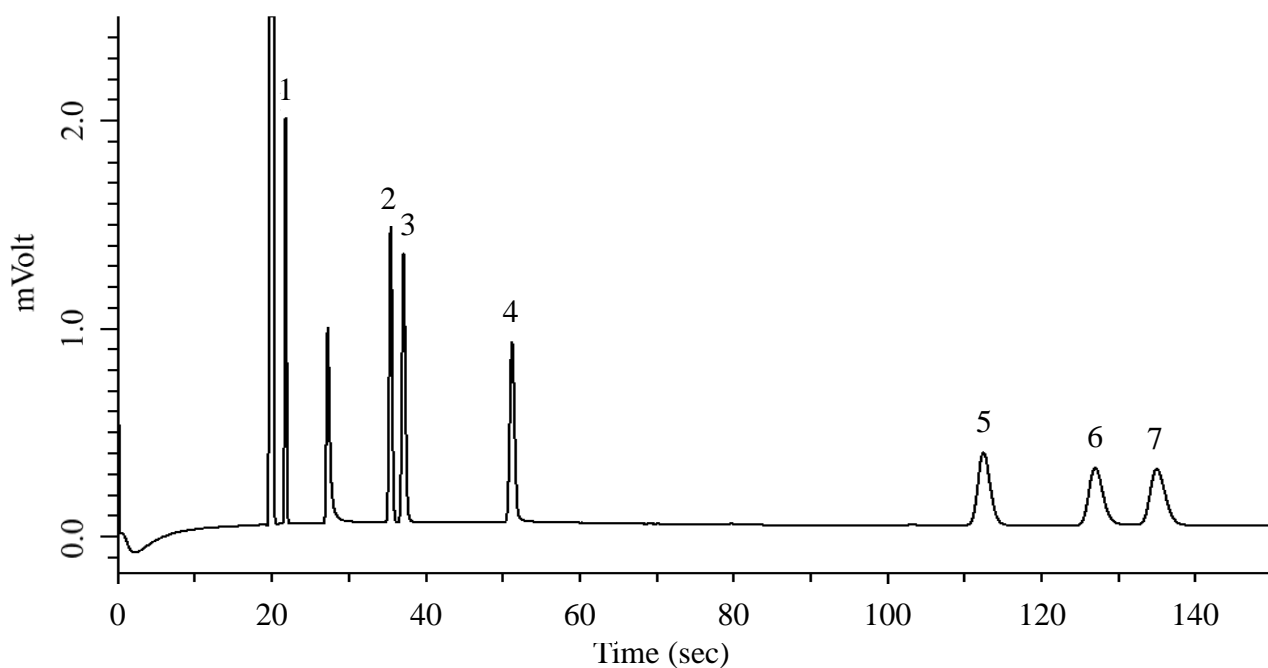

Conditions

System : Micro GC
Column : PoraPLOT Q 10 m
Column Temp. : 120 °C
Carrier Gas : He 190 kPa
Inject Time : 40 msec

Analyte:

1. Carbon dioxide (CO₂)
2. Propylene (C₃H₆)
3. Propane (C₃H₈)
4. Acetaldehyde
5. Propylene oxide
6. Propionaldehyde
7. Acetone