

InertSearch for GC


Micro GC Applications

Inorganic gases (O₂, Methane, CO)

Data No. MA041-1110

Specification

One GC Channel

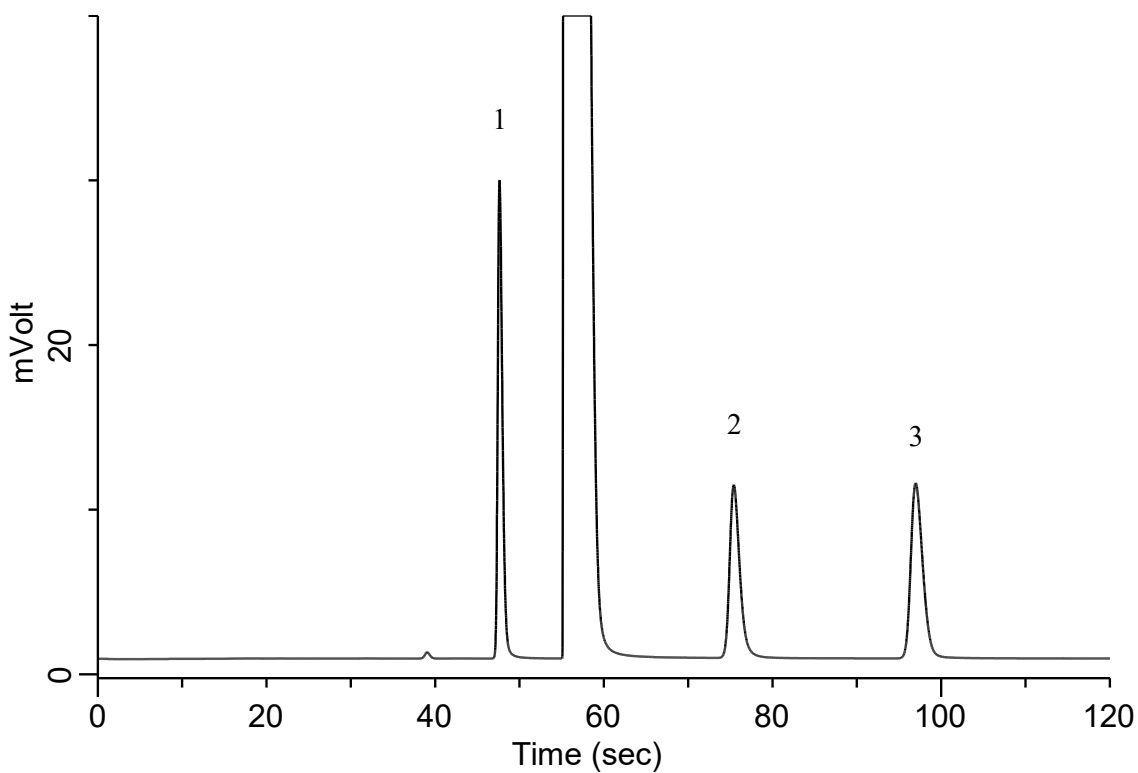

Conditions

System : Micro GC
Column : Molsieve 5A 10m
Column Temp. : 100 °C
Carrier Gas : He 170 kPa
Inject Time : 40 msec

Analyte

1.Oxygen (O₂) 2.5 %(v/v)
2.Methane (CH₄) 2.5 %(v/v)
3.Carbon monoxide (CO) 2.5 %(v/v)

Balance Gas N₂