


Hydrocarbons C3

Specification

One GC Channel

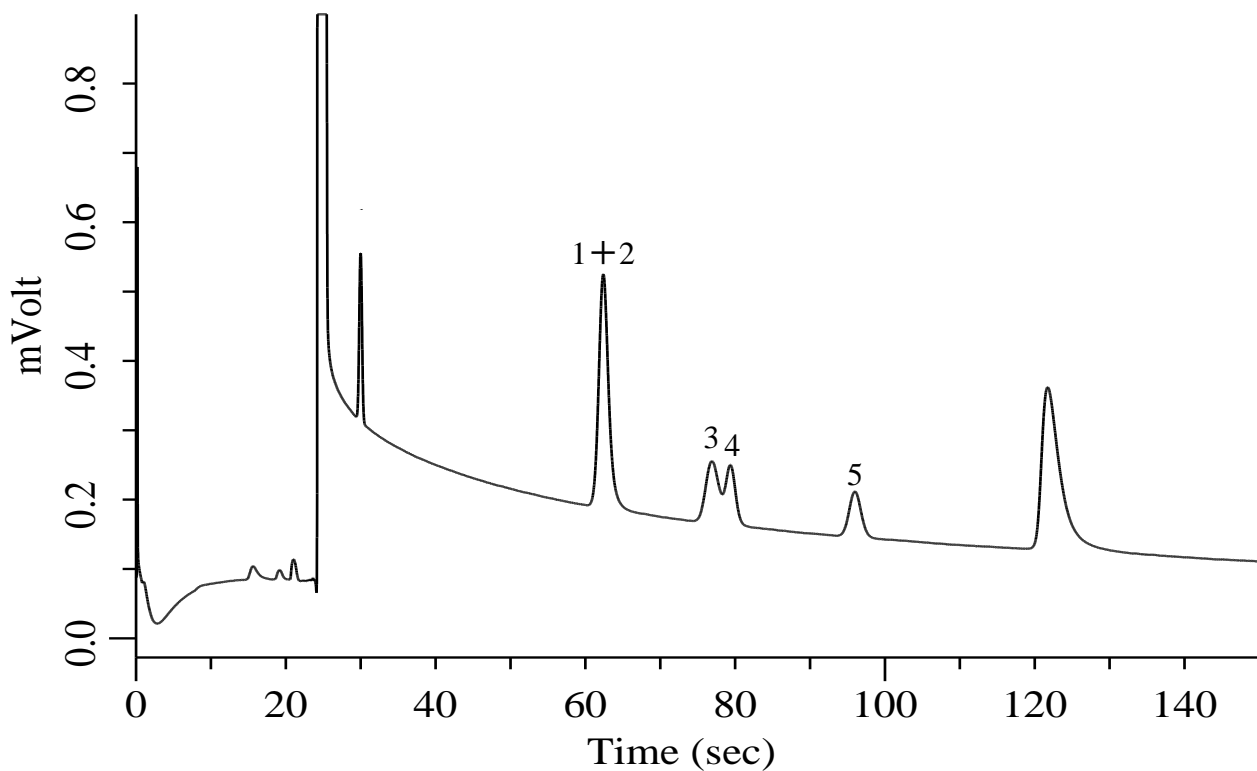

Conditions

System : Micro GC
 Column : PoraPLOT U 10 m
 Column Temp. : 80 °C
 Carrier Gas : He 200 kPa
 Inject Time : 40 msec

Analyte:

1. Propane (C₃H₈) 100 ppm (v/v)
 2. Propylene (C₃H₆) 100 ppm (v/v)
 3. Cyclopropane (C₃H₆) 100 ppm (v/v)
 4. Propadiene (C₃H₄) 100 ppm (v/v)
 5. Methylacetylene (C₃H₄) 100 ppm (v/v)

Balance Gas N₂